

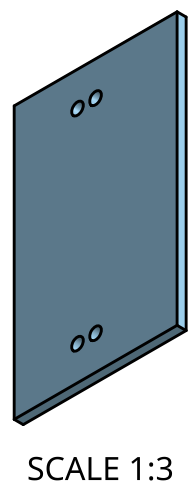


B

A

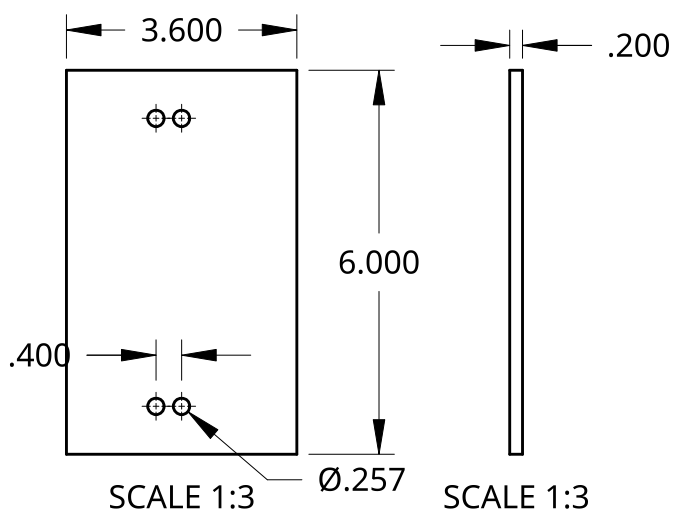

The backplate is tightened to the face of the axles to stop rotation.

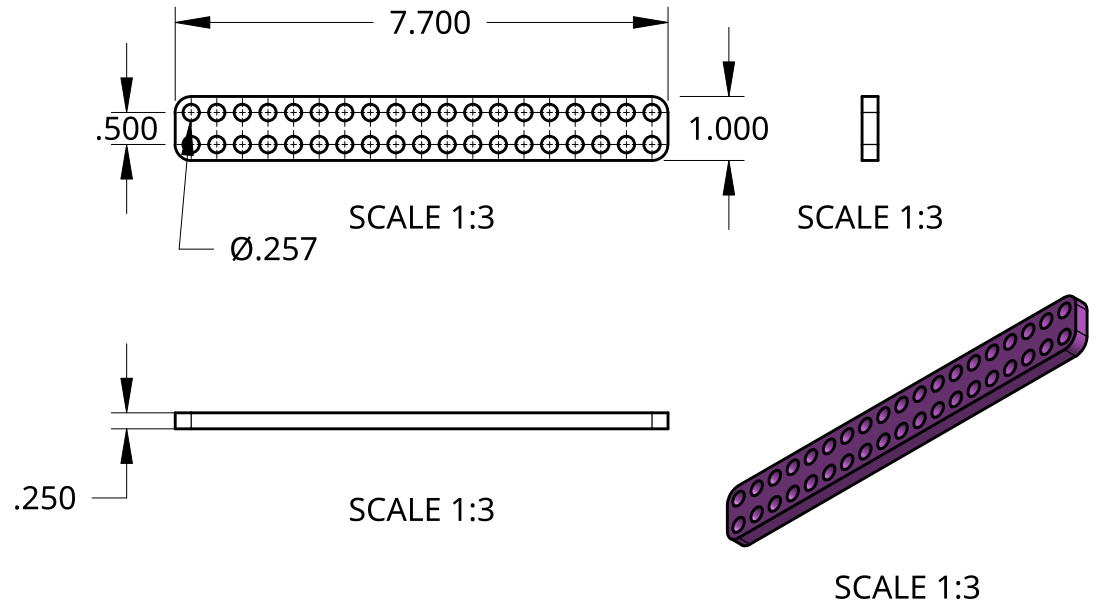

The extender is used to move the servo closer or further to the ground relative to the elevation of the axles' intersection.

These pieces are used in conjunction to the main card style piece to better attach and secure the servo to the bottom/back of the wheelchair. The holes are compatible with 1/4 inch screws and are recommended to be tightened with wingnuts for easy installation.

UNLESS OTHERWISE SPECIFIED, DIMENSIONS ARE IN INCHES . XX = ±.0- .XXX = ±.00- .XXXX = ±.000- SURFACE FINISH  DO NOT SCALE DRAWING BREAK ALL SHARP EDGES AND REMOVE BURRS THIRD ANGLE PROJECTION 		NAME	DATE	TITLE Card Style Solution: Accessory Pieces Mockup 1		
	DRAWN	ROLI VELA	04/27/2026			
	CHECKED	CEO of IDEO: Michael Peng	4/30/2026			
	APPROVED	Prof. Lake	5/1/2026			
				SIZE A	DWG NO. 1	REV. 1
	MATERIAL PLA	FINISH n/a		SCALE 1:2	WEIGHT Var.	SHEET 1 of 2

2

1